

# HemoCue<sup>®</sup> Hb 801 System



## Innovation at its best from the pioneers in Hb testing

### Hb 801

With the introduction of HemoCue<sup>®</sup> Hb 801 System we are taking the Hb point-of-care testing to a new level. The system is designed with your needs in mind, maintaining the robustness and quality you expect from HemoCue products. It provides ease of use and enables work flow efficiency with connectivity options delivering fast result at the time of care.

HemoCue<sup>®</sup> Hb 801 System is optimized for hemoglobin testing in hospital, primary care, public health and road ambulance settings.\* The system provides quick, easy access to lab-quality results without compromising accuracy, even in demanding climates with high temperatures and humidity.

### Maximize time with your patients

- Immediate result in less than 1 second
- Work everywhere with flexible power supply options
- Easy to clean - less down time

### Experience daily operational efficiency

- Large display and symbols allow you to quickly monitor workflow and result
- "Let-go-grip" - no need to keep holding the microcuvette down during measurement
- Vendor-neutral connectivity offers flexibility for e.g. data transfer

### Trust your results at the point of care

- Factory calibrated against the ICSH recommended reference method. No further calibration needed
- Robust design and accurate testing within a wide range of temperatures and humidity
- Built in self-test

*\* Hospital and road ambulance settings only for markets where EU IVDD CE mark is applicable*

# HemoCue® Hb 801 System

**Specifications**

**Principle**

The HemoCue® Hb 801 System consists of an analyzer together with microcuvettes. Blood is drawn into the microcuvette cavity by capillary action. The measurement takes place in the analyzer, which measures the absorbance of whole blood at an Hb/HbO<sub>2</sub> isosbestic point.

**Calibration**

The HemoCue® Hb 801 System is calibrated against the hemiglobincyanide (HiCN) method, the international reference method recommended by ICSH, for the determination of the hemoglobin concentration in blood. The system is factory calibrated and needs no further calibration.

**Sample material**

Capillary or venous whole blood

**Measurement range**

1.0-25.6 g/dL (10-256 g/L, 0.62-15.9 mmol/L)

**Measuring time**

< 1 sec

**Sample volume**

10 µL

**Dimensions**

143×87×45 mm (5.6×3.4×1.8 inches)

**Weight**

< 250 g (0.55 lbs), batteries excluded

**Storage temperature**

Analyzer: 0-50 °C (32-122 °F) Microcuvettes: 10-40 °C (50-104 °F). Opened or unopened vials of microcuvettes can be stored for a shorter period of time (6 weeks) between -18-50 °C (0-122 °F).

**Operating temperature**

10-40 °C (50-104 °F)

**Power**

USB adapter or 3 AA batteries (1.5 V Alkaline or 1.2 V NiMH).


Optional: HemoCue Rechargeable battery

**Interface**


USB and Bluetooth Low Energy  
Vendor-neutral connectivity

**Quality control**


Built-in self-test, optional liquid controls

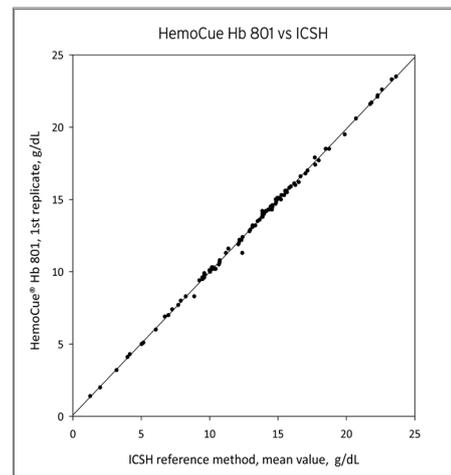

## Three simple steps

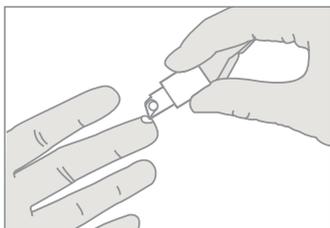

**1** Fill microcuvette.

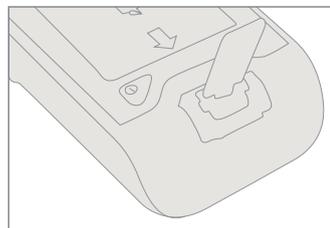

**2** Place microcuvette into analyzer.



**3** View results.